


CENTRELLA[®] SMART+ BED TIPS




GRAPHICAL CAREGIVER INTERFACE (GCI)[®] TOUCHSCREEN

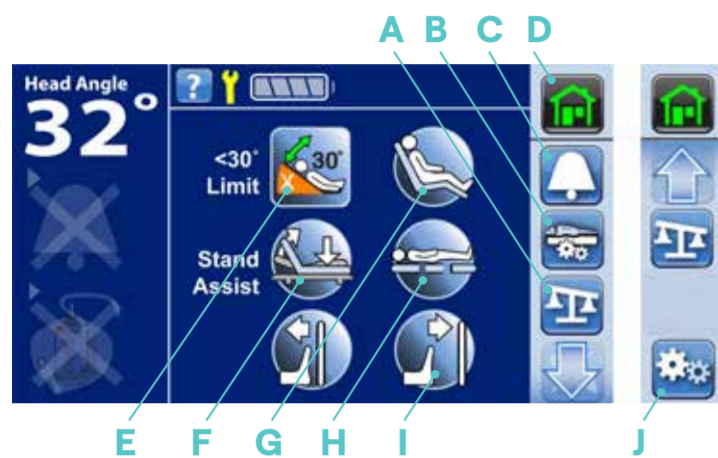

Controls	
A	Scale
B	Air Mattress Controls*
C	Bed Exit Alarm
D	Home Button
E	<30° Head-of-Bed Angle Limit
F	Stand Assist
G	Dining Chair [®] Function
H	Bed Flat Position
I	FlexAfoot [™] Feature
J	Settings/Preferences

Weighing the Patient:

- Select the Scale button on the Navigation Menu on the GCI touchscreen.
- Select Weigh and follow the on screen instructions then select Continue.
- Select Save to save the patient's weight. You will hear a single beep confirming the weight has been saved.

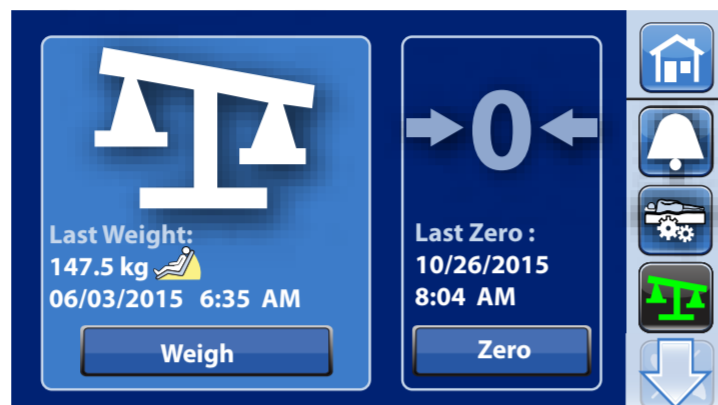

STAND ASSIST

To Activate Stand Assist:

- Select the Stand Assist icon on the home screen of the GCI touchscreen.
- Select and hold the Stand Assist control until the legs are flat and the head is at the desired angle.
- One-button simultaneously levels legs to a flat position and raises the head-of-bed to position the patient safely.

Note:

- If equipped with a Centrella[®] max powered air surface, the seat section of the bed will automatically inflate to provide more support for the patient.

Pro Tip:

- Teach your patients to use Stand Assist on their patient pendant*. Non-fall risk patients can use this on their own and fall-risk patients can use this proactively before you get bedside. *Patient pendants are an optional feature.

FOLEY POSITION LIMIT



- The Foley Position Bag Holders can be found at the foot-end of each side of the bed.

To Turn the Foley Position Limit On and Off:

- Select the Setting Menu control on the Navigation Menu of the home screen and select the Bed Features control.
- Select Off or ON, and Accept to lock in the new setting.

Note:

- A Foley Position Limit indicator will display on the screen, when activated.

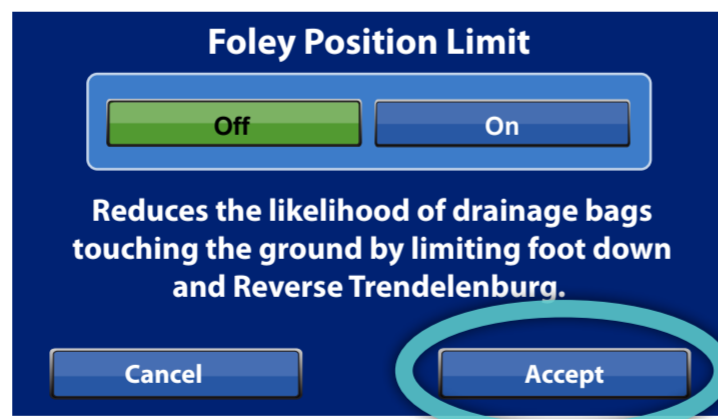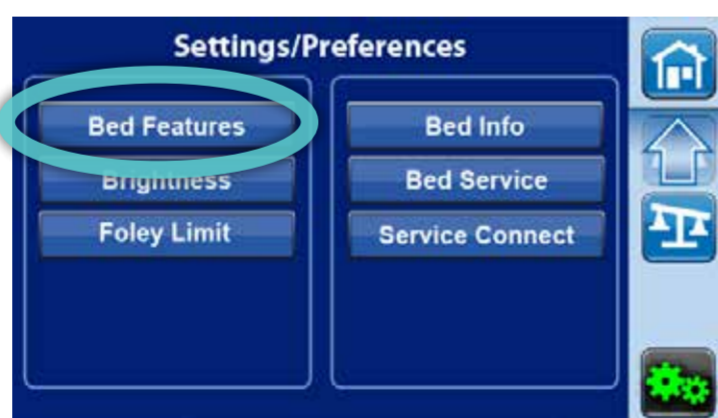

CENTRELLA BED EXIT ALERT SYSTEM

Three sensitivity settings:

- Changes Position
- Moves Towards the Edge
- Has Left the Bed

To turn on the bed exit alert:

- Select the Bell icon in the Navigation Menu of the GCI touchscreen.
- Select the appropriate setting for your patient.
- Listen for a single beep and look for the Bed Exit On icon indicating the system is armed.

Silence the Bed Exit Alert System without deactivating the system.

- Select the silence button on the GCI or yellow silence button on siderail. The alarm will re-arm automatically, once patient is back in bed.



SCALE OPERATIONS

Zeroing the Scale:

- Select the Scale icon on the Navigation Menu of the GCI touchscreen.
- Select Zero. Select "New Patient" or "Re-Zero" for same patient.

New Patient Zero:

- Clears previous patient weight history
- Zeros the scale
- Resets bed to facility's default settings

Re-Zero:

- Keeps patient weight history and bed settings customised for that patient.
- Zeros the scale

AUDIBLE ALERTS

Two Types of Audible Alerts:

- Single Beep:** Activity was performed successfully (i.e., when the bed reaches the lowest position).
- Triple Beep:** Attention is needed; look at the message on the GCI Touchscreen for more information (i.e., when you attempt to use control that is locked-out).


SIDERAIL PROTOCOL STATUS

To Adjust Siderail Protocol:

- Select the Setting Menu control on the Navigation Menu of the home screen and select the Bed Features control.
- Select Rail Protocol and configure rail protocol followed by your facility and select Accept.
- The siderail protocol is reset to facility default when the bed is zeroed for New Patient.

Note:

- The siderail protocol is based on your facility's policy.

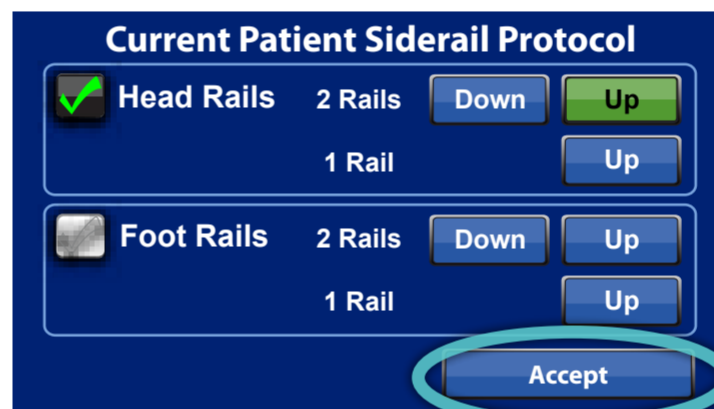

SAFEVIEW[®]+ SYSTEM

Siderail Protocol:

- Green = Number of siderails up is equal to, or more than, set protocol
- Flashing Amber = Fewer siderails up than set protocol

Bed In Lowest Position Status:

- Green = In lowest position
- Flashing Amber = NOT in lowest position

Bed Exit Status:

- Green = Bed Exit is ON
- Flashing Amber = Bed Exit ALERTING
- Solid Amber = Bed Exit SILENCED
- Blue = Bed Exit is OFF



INTELLIDRIVE[®] POWERED TRANSPORT

To Activate IntelliDrive[®] Powered Transport:

- Unplug the power cord and set the brake pedal to Steer (green pedal down).
- As the drive wheel lowers, the indicator on the Transport Pod flashes green. When the wheel has lowered completely, the indicator turns solid green and a single beep will sound.
- Select down on one of the black enable switches on the underside of the blue push handles. Apply pressure to the handles to move the bed either forward or backward.

To Deactivate IntelliDrive[®] Powered Transport:

- Set to Brake or Neutral, or plug the bed into power.

AIR MATTRESS CONTROLS*

Three Mattress Modes:

- Normal Mode:** Provides continuous full-body pressure redistribution for patients.
- Max Inflate Mode:** Maximises firmness of the surface, for 15 minutes, to assist in patient handling and repositioning.
- Turn Assist:** Helps the caregiver turn the patient to the left or right. This can be used for procedures such as patient positioning and linen changes.

Navigate to Different Mattress Modes:

- Select the patient mattress icon.
- Select the most appropriate setting for your patient.

Note:

- For proper operation, there must be a minimum of 32 kgs (70 lbs) on the surface.

Pro Tip:

- Press and hold Adjust Bed Button to adjust the head-of-bed, knees and foot angle to the required position.



Hillrom[™]

To learn more call Hillrom Customer Care on 1800 445 576 or visit hillrom.com.au

*Optional Feature

Learn more at hillrom.com.au

Always read the labels and instructions for use. Contact your Hillrom Representative for the latest version.

AUTHORISED AUSTRALIAN SPONSOR Hill-Rom Pty Ltd. Suite 4.01, 2-4 Lyonpark Road, Macquarie Park NSW Australia 2113
© 2021 Hill-Rom Services, Inc. ALL RIGHTS RESERVED. APR206701_AU-r1 03/2021

FOR DETAILED WARNING AND CAUTION STATEMENTS, AND OPERATING INSTRUCTIONS, REFER TO THE CENTRELLA[®] BED USER MANUAL (193587).